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
# Instruction

- Can be worn in the lower or upper jaw
- The integrated wire should be bend easily according to the dental arch – both water pads should remain passively and symmetrically between the teeth
- The V-mark shows the midpoint of the connection tube of both water pads. It should be rested/located near the midline of the upper or lower jaw
- Please wash up the **AquaSplint®mini** with water after removal and place it into the hygiene box

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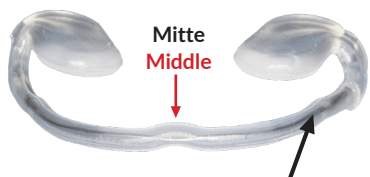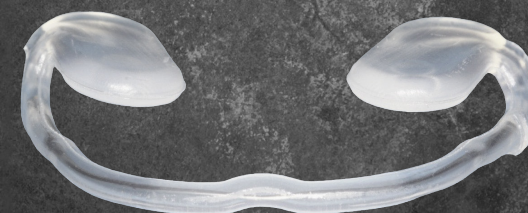
### Application/Indication:

**AquaSplint®mini** is a self-adjusting bite-splint that is filled with water and offers immediate therapy to the TMD syndrome and associated disorders like TMJ-pains, migraine, cervical syndrome, tinnitus, head and face ache, blockades, moderate sleep apnea.

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### Contraindications:

- Severe bruxism/grinding reduces the durability
- Due to sharp edges of prostheses, dental crowns etc. the water pads could be damaged
- Missing molars - at least one second premolar tooth in every quadrant should be present

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### Characteristics of AquaSplint®mini:

- Lower height – thus particularly applicable for TMD-patients with open bite
- Preferred by non-dental physicians like orthopedists, ENT physicians, also osteopaths, physiotherapists and pharmacists
- No relining
- Allows balancing of the lower jaw
- Neutralized any forced bite
- The integrated wire can be bent individually pursuant to the dental arch
- Effective and immediate access to therapy
- With progress and pain relief an occlusal rehabilitation (dentist) or an **AquaSplint®thermo** can follow if needed

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### Advantages compared to classic bite splints:

- Quick pain alleviation and muscle relief
- Installable instantly without lab work
- Universal Size
- No sharp edges
- No more than 10 hours of daily wearing
- No appointments for adjustment (thanks to the hydrostatic balance between both water pads)
- Alternative to both relaxation and distraction splints.
- Can be used during orthodontic treatment

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### Difference between AquaSplint®thermo & mini:

**AquaSplint®mini** is about 40% thinner than the **AquaSplint®thermo** due to the omission of both saddles and relining. Therefore it is more suitable to treat open bite patients. **AquaSplint®mini** has in comparison to the **AquaSplint®thermo** an integrated metal arch wire in the tube which is connecting the two pads. This wire makes it possible to be individually adapted to the shape of the dental arch.