

TROMATZ®  
Simple Pro  
Rechargeable

### How to use



1

Apply toothpaste



2


Press button to activate  
microcurrent waves




3

Brush as you like  
until the blue light  
turns off

### How to replace brush heads

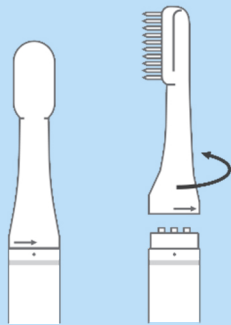


Turn counterclockwise  
remove used brush head

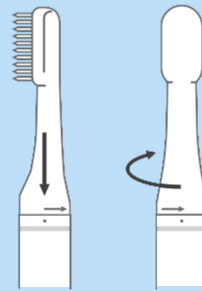

Align new brush head  
and turn clockwise

\*We recommend replacing brush heads every 3-4 months

### How to charge

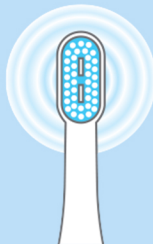

Charging required

LED will blink  
when battery  
is less than 20%




Charging

Connect C-type  
cable to toothbrush  
**Charging port  
must be completely  
dry before charging**



Charging

Blinking LED will  
indicate charging  
It takes 3 hours to  
fully charge.

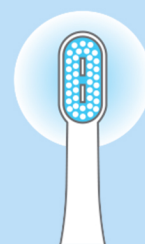

Charging complete

Fully charged LED  
will stop blinking  
And it can last  
up to 2 weeks.