Innovative Idea and Possibility **Good Doctors**


Intraoral Camera





Good Doctors Co., Ltd.

GoodDrs Camera Classification

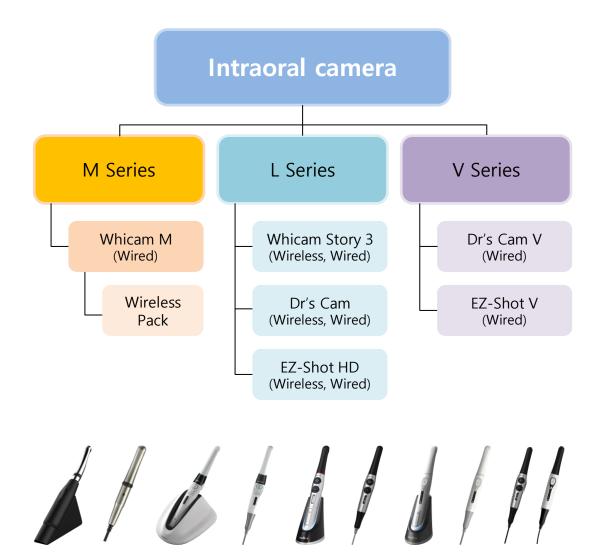



(Multi Functional station) All cameras can be used with station

Strong point of GoodDrs Camera



Auto focus Real Autofocus function shows clear images.



HD image High definition show more detail images,



Software compatibility

Can be use camera with most of dental imaging software, * It would need twain driver for some of dental imaging software,





Zoom function

Ability to enlarge the image with zoom feature to support user convenience.



* Need to install "Zoom function" program on the PC.






Good Doctors Co., Ltd.

Strong point of GoodDrs Camera

Multi Twain Driver

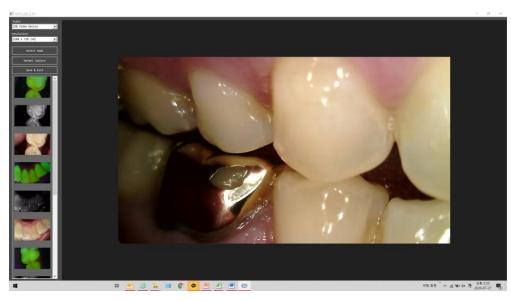


HID Manager

HIDDrsManager Ver2.0.2	- • ×
KEY-VALUE Space Enter Space->Enter Space->Space Enter->Space Space+S->Space F5->F5 Ctrl+L T->Space G->Enter Alt+F->Alt+L Alt+F->Alt+G Enter->LCtrl+LShift+V Command(WIN)+[Mouse L Mouse D Advanced key delay Current delay(ms): 500	Device Device connection status Current Custom-Key SPACE Custom Capture-Key Custom KEY-VALUE Apply Close O Apply
Heat Time	
Time(min): Off Off	On

∕]₽_{,,,,,,},

Remote Control

원격지원 서비스 Remote Assistance Service

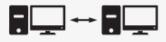

Chanel Chanel 2



Tel: +82 32 424 6325 E-mail: info@gooddrs.com

Features of Wireless camera

Wireless type only

(((i)))

High speed digital wireless High speed digital wireless provide more clear live image without interruption.



Wireless receiver Compact wireless receiver without any power supply.



Big capacity battery

Long time use with big capacity battery,







Whicam M

Whicam Story 3 CS

Dr's Cam

EZ-Shot HD

Receiver





M-Series

Macro to Infinity Intraoral Camera Model Name : Whicam M

- Durable metal body
- It is the only product on the market that can obtain Macro to infinity Image
- · Compact size with a dental handpiece design which is familiar to dentists
- · Can be easily mounted on a dental unit chair
- · Various color impression functions and convenience functions such as a built-in mouse

Whica	am Macro to Autofocus Dente	O Infinity al Cemera with Mouse		$\overline{\mathbf{n}}$
Macro to in	hfinity	Dental handpiece design		
AF Auto-focus	s Twain	Support twain driver		
METAL Metal body	/ case <u> t</u> <u> t</u>	Slim head size		
Mouse fund	ction FHD	FHD image support		
			Full HD 1080p	Whicam M

M Series

Whicam M



Zoom function Ability to enlarge the image with zoom feature to support user convenience.

LED brightness adjust Ability to adjust the LED brightness.

Support twain driver Support twain driver allows for use in almost any dental imaging software.

Software compatibility Able to be used with almost any dental imaging software in the market.

Detachable USB cable Detachable USB cable makes it easy to move the camera to other dental chairs.

=









FHD image support FHD The device shows FHD (1080p) in live image supported by a 5M sensor.








Color impressions Adjust the color impression on camera to fit different monitors or screen.







Same color chart image with switte

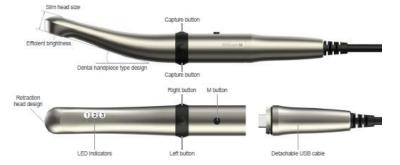

Wireless Pack



- The wireless pack make the product to be wireless Intraoral camera
- It can be supplied as option for existing wired type, or can be supplied wireless type
- The wireless technology applied to products is proven technology used over many years


M Series



L-Series

Wireless and Wired Intraoral Camera System and Extraoral Illumination Optic System

- · Possession of a variety of models of wired wireless intraoral cameras
- Auto focus on all products
- Wireless products are the first digital wireless in the market (High speed digital wireless without noise)
- User-friendly features such as a mouse
- Own patient image management software

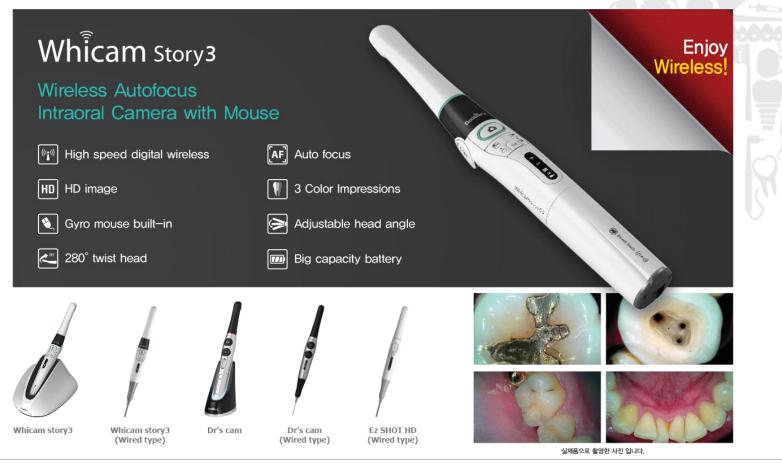

Whicam story 3



Wireless type only



16

High speed digital wireless High speed digital wireless provide more clear live image without interruption.

Wireless receiver Compact wireless receiver without any power supply.



Big capacity battery Long time use with big capacity battery,

Auto focus AF Real Autofocus function shows clear images.

HD

HD image

High definition show more detail images,

Adjustable head angle User can use the camera with prefer angle,



Gyro mouse built-in Useable as a mouse,



280° twist head Easy to get approach anywhere in oral cavity,



OLED screen

Bright and clear OLED screen indicates various device status,



Adjustable LED brightness Can adjust LED brightness to avoid light reflecting,



3 Color Impressions Can change live image color impression (Original, Warm, Cool)



Software compatibility

Support keyboard short-cut data changing program and twain driver,

Dr's cam



Wireless type only

High speed digital wireless High speed digital wireless provide more clear live image without interruption,

Wireless receiver

Auto focus

Can get clear image with auto focus function in oral cavity,



(((*)))

AF


280

OLED

淤

Can get high definition image without distortion.

Gyro mouse built-in Useable as a mouse.

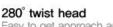

Easy to get approach anywhere in oral cavity,

OLED screen Bright and clear OLED screen indicates various device status,

Adjustable LED brightness To avoid image quality loss by light reflection, can adjust LED brightness,



Software compatibility

3 Color Impressions

Can be use camera with most of dental imaging software, * It would need twain driver for some of dental imaging software,

Can change live image color impression (Original, Warm, Cool)









EasyShot HD



Wireless type only



High speed digital wireless High speed digital wireless provide more clear live image without interruption,



Auto focus

Can get clear image with auto focus function in oral cavity,



AF

HD image Can get high definition image without distortion.



Adjustable LED brightness To avoid image quality loss by light reflection, can adjust LED brightness,



Black & White function Can take and send analog X-ray film image to PC.

Plug & Play

No need additional driver to install the device,



Easy to capture

Easy to image capture by camera only without extra device like foot pedal,



Software compatibility

Can be use camera with most of dental imaging software, * It would need twain driver for some of dental imaging software,




V-Series

Current Product

V series Intraoral Camera Model Name : Drs Cam V, Ez SHOT V

- Product that increases price competitiveness
- Competitive price and high definition
- Convenience functions such as mouse

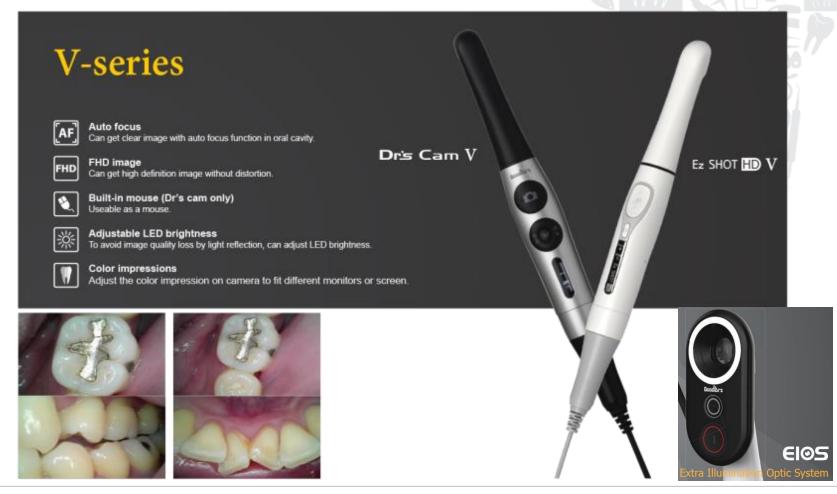

V Series

M station

Multifunctional Station for Dental Camera All in One System without Computer Model Name : M station

- Competitive price product
- Targeting Normal Dental office and Unit chair companies
- Management and consultation of most patient images without a computer
- Easy installation



M station



Features of M-STATION



Without computer system


It has its own image management program, Patient image management is possible without PC or S/W



Snapshot & Video Capture Snapshot & Video Recording is possible



High definition (HD) image support High-resolution display supported via HDMI port

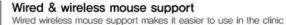



Various image editing functions

Various editing functions such as zoom in, zoom out, mirroring, drawing for communication with patient








USB memory card Copy saved data using USB memory card



Various language support Supports English, German, Korean, Chinese, Italian, Spanish

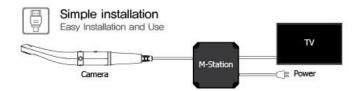

-



Supports three UI themes

Selectable UI style for user's desire












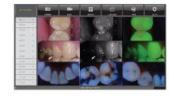





Download 1,000 possible files Store up to 1000 files on internal memory



Compact size 82 x 82 x 17mm / 96g

Extra Intraoral Optic System (EIOS)

- · Can take full arch or face with EIOS
- Can take external images using strong lighting
- · Easily detachable with all cameras of GoodDrs

EIOS

Extra Ilumination Optic System



Extra oral image Help to show full arch photo with Good Doctors intraoral cameras



Compact design Compact & Easy to apply design



Rechargeable battery User doesn't need to change battery



Bright LED 24 LED make light



NEV ARRIVA



Current Product

Accesories



Disposable Cover

35*100(S) - Smaller 35*100(S) - Wider 46*195(S) - Wider 60*205(S) - Wider



AP.

Ľ١



Receiver



Cradle



Comparison Table



Comparison table with Good Doctors intraoral camera

Model	Wireless / Using time	Focus Type	Focusing Range	Image Quality (Max)	Mouse	Casing	Capture Button
Whicam M	Yes / 2 hours	Auto Focus	5mm ~ Infinity	FHD (1980 x 1024)	Yes	Metal	Dual
Whicam Story3	Yes / 2 hours	Auto Focus	5mm ~ 50 mm (full arch)	HD (1280 x 720)	Yes	Plastic	Single (Rotate head)
Dr's Cam	Yes / 25 mins.	Auto Focus	5mm ~ 50 mm (full arch)	HD (1280 x 720)	Yes	Plastic	Single (Rotate head)
EZ Shot HD	Yes / 25 mins.	Auto Focus	5mm ~ 50 mm (full arch)	HD (1280 x 720)	No	Plastic	Single (Rotate head)
Dr's Cam V	No	Auto Focus	10mm ~ 15 mm	FHD (1980 x 1024)	Yes	Plastic	Single (Rotate head)
EZ-Shot V	No	Auto Focus	10mm ~ 15 mm	FHD (1980 x 1024)	No	Plastic	Single (Rotate head)

* Whicam Story3, Dr's Cam EZ Shot HD the focusing length is extended to infinity in case EIOS apply

* Dr's Cam V, EZ Shot V the focusing length is extended to full arch in case EIOS apply





Thank You



Good Doctors Co., Ltd Head Office :#208, B-dong, 283 Bupyeong-daero, Bupyeong-gu(Woolim Lions valley, Cheongchoen-dong) Incheon, 21315, Republic of Korea