



# Dentozone CATALOG



NEW COLLECTION



## GET IN TOUCH!

### Location

#1011, Ace High-end Tower  
1 cha, 5, Digital-ro 26-gil,  
Guro-gu, Seoul, 08389 S.KOREA

### Web

WWW.DENTOZONE.COM

### Phone Number

+82-2-868-2875

### E-mail

MARKETING@DENTOZONE.COM

WWW.DENTOZONE.COM



# Welcome to product Catalog

Dentozone is the first company in the world to have developed the sterilizing dental water system, a dental equipment of the next generation, and is always striving to develop innovative products and satisfy customers.

Dentozone's sterile dental water system developed and manufactured by Dentozone Corporation has proven safety, effectiveness and convenience for dental professionals that chemical disinfectants or similar products from other companies cannot match.

Moreover, it is an eco-friendly product that does not harm humans or the environment and can most effectively purify and sterilize dental unit water that are quite costly.

The stylish and modern design, various high-tech functions available for customers' convenience, and the improved customer service all make our products stand out in the international dental market. Also our products have been certified by many well-recognized international associations, and can guarantee quality and safety.

Dentozone will continue to research and develop in order to achieve our customer's satisfaction with a wide range of choice.

We ask for your continued support while Dentozone further grows into the leading company in the world's dental equipment market.


# Table Content

	<p><b>eBite3S Blue light</b></p> <p>eBite3S Blue light can check orthodontic attachments easily.</p>	<p><b>4-5 PAGE</b></p>
	<p><b>eBite3S</b></p> <p>eBite3S is a rechargeable cordless loupe light.</p>	<p><b>6 PAGE</b></p>
	<p><b>Portable eBite</b></p> <p>As a cordless type intraoral illuminator, you can see your mouth brightly without shadows.</p>	<p><b>7 PAGE</b></p>
	<p><b>99 Suction Block (Reusable)</b></p> <p>It is a 99 suction block that automatically suctions without the assistance of an assistant during treatment.</p>	<p><b>8-9 PAGE</b></p>
	<p><b>99 Suction Block (Disposable)</b></p> <p>This is a disposable type. Suction is performed automatically without the assistance of an assistant.</p>	<p><b>10 PAGE</b></p>
	<p><b>DPS2</b></p> <p>Dentozone Corp, the No. 1 dental water quality management company in Korea. It manages dental hospitals nationwide in Korea and manages dental water quality with plasma.</p>	<p><b>11 PAGE</b></p>

# New orthodontic Blue light

The eBite3S blue light distinguishes the resin, which was difficult to distinguish with the naked eye during removing resin, as a fluorescent blue light, so it is a necessary light for delicate work that can preserve natural teeth without erasing them during resin removal.



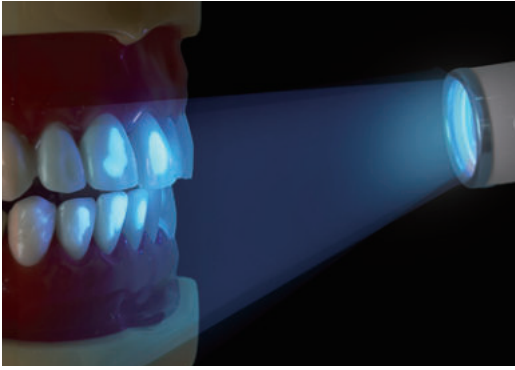

## Distinguish transparent orthodontic attachments

When removing orthodontic attachments that were difficult to distinguish due to their similarity to the color of teeth are distinguished by a **fluorescent blue color** and can be easily removed.



## Attaches to goggles, loupes, and glasses

It can be clipped to goggles, loupes, and glasses, so both hands are free when removing the attachment.



## Angle can be adjusted up to 90°

The angle can be adjusted down to 90°, so you can see the treatment comfortably from any angle.



## Up to 5 hours of use

It is a rechargeable type product that can be used for up to 5 hours on a single charge.



### • Product specification

Product name	eBite3S Blue Light
Size	41mm X 70mm X 54mm
Weight	41g
Brightness step	Step 3
Brightness	(Low) 1,300Lx (Mid) 1,700Lx (High) 2,000Lx
Charging time	2Hours
continuous use available time	Continuous use in weak mode: up to 5 hours Continuous use in strong mode: up to 3 hours


# eBite3S

Rechargeable  
cordless  
LED illuminator



**Fills the entire field of view from edge to edge.**


Featuring eBite3S illumination with accurate color rendering that fills the entire field of view from edge to edge.

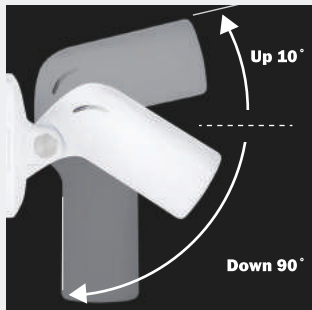

**Adjusting the angle**



**Attach the goggles**

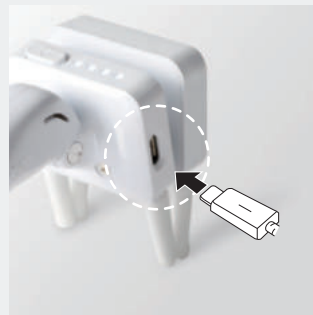

**Charging the Battery**



**Curing Filter**



• **Product specification**

Product name	eBite3S
Input	AC 100~240V / 50-60Hz
Output	DC5V · 1A
Charging time	2Hours
Weight	39.9g
Illumination	High 6,000Lx / Mid 4,800Lx / Low 3,000Lx
Color Temperature	6,500K
Life Span	50,000 Time



portable  
**eBite**

Rechargeable cordless  
Intraoral Light



**Illuminates the inside  
of mouth without shadow**

Cordless & Rechargeable illumination solution for the oral examination

**Cordless & rechargeable illumination**  
solution for oral examination

**LED light close to natural light**  
Life of LED is 50,000 hours

**Designed for optimal patient comfort**

**Designed for optimal patient comfort**

**Bite block and LED cover can be disinfected through autoclave.**  
Normally working temperature 200°C  
(Disinfect them for 3 minutes in autoclave of 135°C)

**Combining bite block with eBite It spreads the mouth of patient naturally**



• Product specification

Product name	Portable eBite
Input	AC 100~240V / 50-60Hz
Output	DC5V · 1A
Charging time	1Hour (60min)
Illumination	5,000~8,000 Lux
Color Temperature	5,500K
Life Span	50,000 Time

DENTOZONE

## 99 Suction Block®


Reusable



# New Self Suction 99 Suction block

The ergonomically designed 99 suction block gently wraps around the gums and removes more than 99% of water and saliva from the mouth during treatment in real-time. It also removes saliva inside the throat, which was difficult to remove even with conventional suction tips, making patients and patients comfortable.

**Installation** is **super easy!**

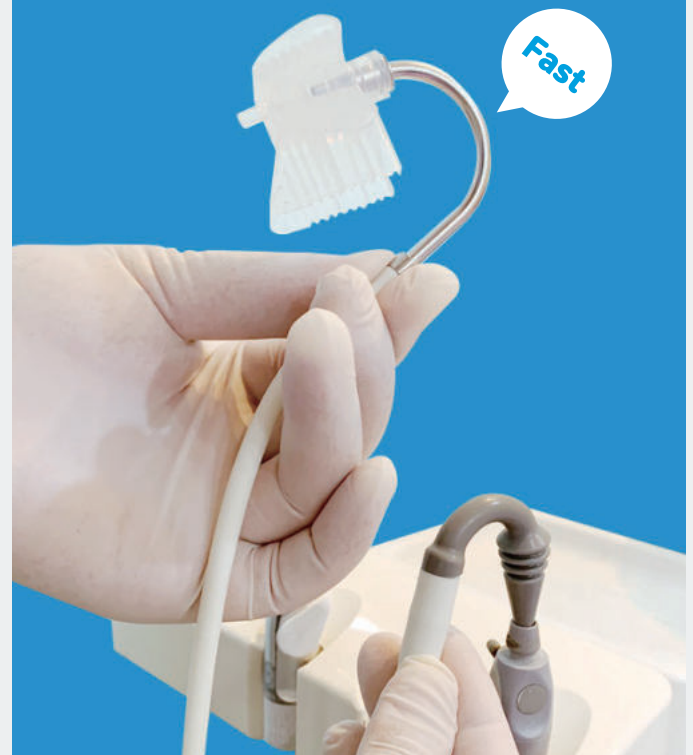



## 99 Suction Block

It was developed with the convenience of both patients and staff in mind.



It's not hard, it's soft.

### 99 Suction Block made possible for 99% suction!

99Suction Block sucks water out of patient's mouth continuously as it wraps gum of teeth softly.

### Made to work for both mouth gag and suction.

It does not tire patient out for a long treatment. Its suction capability works seamlessly without stop for treatment and assistant needs not interrupt the treatment with the use of suction tip.

### Reusable

99 suction blocks can be reused after sterilization with autoclave. It's very convenient because you can use it again after using it.

### Not hard and soft.

The 99 suction block made of soft silicon is not hard when the patient bites the product, so there is no burden on the oral cavity.




#### • Product specification

Product name	99 Suction Block(Reusable)
1 Package	Suction Block Connector(1ea) / Suction Block (5pcs)
Package size(mm)	110(W) x 100(L) x 100(H)
Weight(g)	177g
1pcs(size · weight)	32(W) x 15(L) x 40(H) · 10g

# DENTOZONE 99 Suction Block®

Disposable

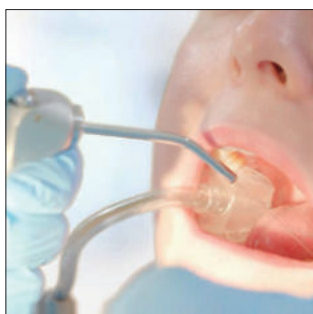

### No need assistant

99 Suction block can suction 99% of the water and saliva in a patient's mouth without the assistant's help.



### Super easy to install

Just Install it into SE.



### Individual Package

It is hygienic because it is individually wrapped.



### Easy to use

There are 20ea in a box, so you just pull out at treatment.




### • Product specification

Product name	99 Suction Block(Disposable)
1 Package	20ea
Package size(mm)	110(W) x 190(L) x 140(H)
Weight(g)	350g
1pcs(size · weight)	90(W) x 194(L) x 17(H) · 15g



# DENTOZONE DPS2

NO Chemicals  
No Bacteria  
No Biofilm

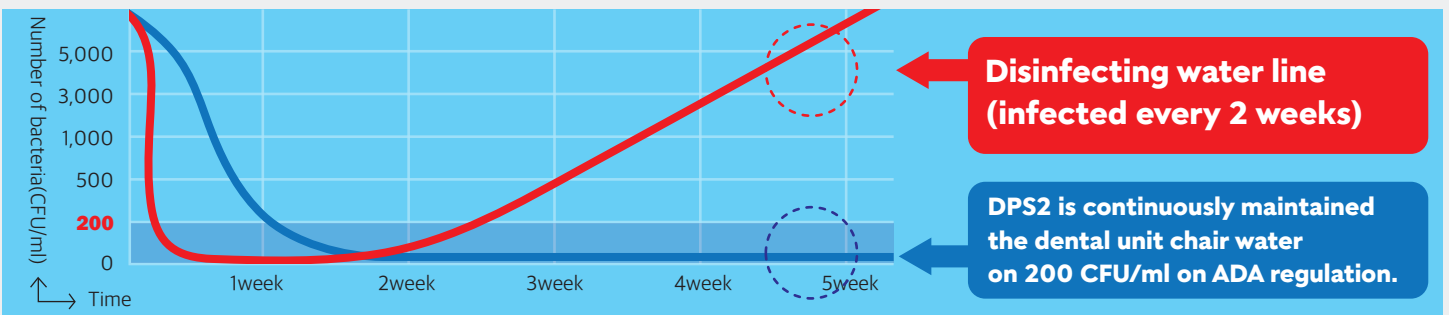

### No chemicals plasma sterilized water

No adding chemicals to plasma water is safe for patients and dentists.

### Interaction with patient on display

The display shows the patient that the water is safe.

### Real time supply of sterilized water

Plasma water supply in real time for using even a drop of water with the micro sensor.

### Compact size easy to install

Compact size of DPS2 is easy to install.



### • Product specification

Product name	DPS2
1 Package	DPS2, Display, Bracket
Package size(mm)	270(W) x 380(L) x 135(H)
Weight(Kg)	3.1kg
Power	DC 24V (100~240VAC)