





Factory 1 | NO.9, 83 Road, Waryong-ro, Seo-Gu, Daegu, South Korea (zip code: 41754) Factory 2 | NO.41 90 Road, Waryong-ro, Seo-Gu, Daegu, South Korea(zip code: 41759) West Daegu Smart Tech Hub Floor B2 NO. 203 **Tel** | +82-53-351-1772 **Fax** | +82-53-351-1775



# The Best Dental Manufacturer & Export Co.


# dental + insight = articulator / arch

It was intended to express Corident's technical skills and professonalism by using the ideal dental structure, thet is an arc, as amotif in develoment and 'The concept that we make the globe smile' contains the philosophy of the company.

Company Name	
Learship	Dong-Hea Ye(president)
Address	NO.9, 83 Road, Waryong-ro, Seo-Gu, Daegu, South Korea (zip code: 41754)
Contact No.	Tel : +82-53-351-1772 Fax : +82-53-351-1775

## Certification



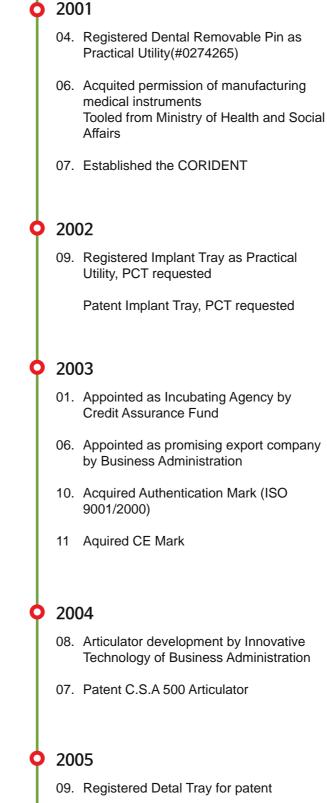

10. Appointed as Venture Company



## 2006

О

Ο

Ο

Ο

Ο

Ο

Ο

П

12. Be selected as Inno-Biz Participation development innovative technology of Business Administration

## 2007

09. Recognized quality management standard for medical devices(GMP)

### 2008

09. Recognition of foreign medical device manufacturers(Japan)

## 2012

03. A dental gold alloy with a tiny amount of titanium patented

### 2014

- 09. A dental occluder patented
- 10. Obtained US FDA accreditation

### 2016

08. Acquired Authentication Mark ISO13485

## 2017

- 01 Cobalt Chromium alloy metal block for dental CAD/CAM and its manufacturing method patented
- 11 The implant Model Articulator patented

### 2019

- 10 Building additional Factory 2
- 11 Recognized quality management standard for medical devices(GMP Class2)

## **CSA 650** ARTICULATOR

THE DIARC650 ARTICULATOR was upgraded from the existing CSA600 Articulator. THE DIARC 650 is an advanced articulator in terms of convenience and design based on ergonomics.



### All adjustable articulator with a calibration system

Condylar Angle(0°~50°) Bennett Angle(0°~45°) I.S.S(Immediate Side Shift)

## Improvement Operation Advancement

The diarc 650 Articulator has veen improved for Condylar Angle, Bennett Angle, Condylar Lock and I.S S(Immediate Side Shift).

Retrusion and Protrution movement are added.

Condylar, Retrusion, Protrution and I.S.S have been advanced equally with human body structure, the design is also refined.

## Convenience






• Fully open accrss from the lingual position for easy access during denture

 Condylar Angle 50° adfusment Latrusion


## | Articulator Development

**Retrusion, Protrusion Movement Reproduction** 

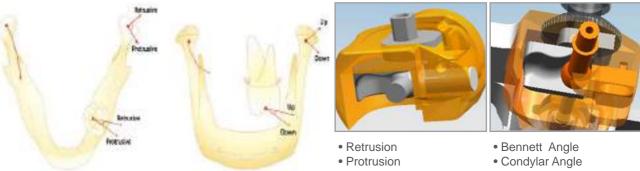

## Improvement of Masticatory Force Control System

Condylar Giudance was improved and reconstrycted with anatomy-based data.

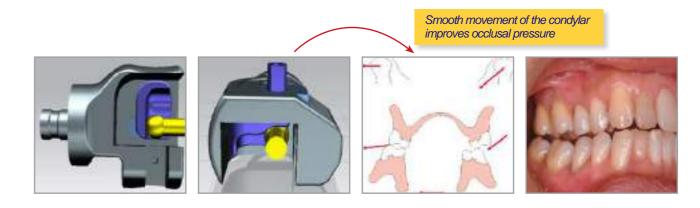





 Retrusion Protrusion



• Simple immediate side shift

## **CSA 600** ARTICULATOR

CALIBRATION System makes it possible for the dental clinics and labs to interchange plaster tooth cast without exchaning articulators themselves. And it ensures that only casts can be iterchanged amoung CSA Articulators without inaccuracy.



### All adjustable articulator with a calibration system

Condylar Angle(0°~50°) Bennett Angle(0°~45°) Side Shift adjustable

## Verstile Articulating System & Higher Precision

System	Calibration System with adhustable condylar inclination, bennett inclination and side shift system.
Precision	CSA600 articulator is the perfect choice for analysis, treatment planning
Operation	It combines the advantages of easy use, simplicity and ability to tilt the instrument back at a 30 degree angle.
Appearance	It is designed to be convenient for work as well as ergonomin, sturdy and lighweight desgin with durable black anodized finish.

## | Convenience





• Convenient attaching & removing cast with Magnetic mounting plate.

• There is no tool required with simply using Accessorty separating screw.

## | Adjustability





- · Adjustable Condylar, Bennett inclination and Side shift system
- INCISAL GUIDANCE : Incisal table 10°&20° and Incisal pin

## | Calibration



Precision, Convenience, **Time-saving** 





• Fully open access from the lingual position for easy access during enture setup so excellent to check lingual position when set in array.



· Handy hard box to carry



• Easy use to tilt the instrument back at a 30 degree angle.





## **CSA 400** ARTICULATOR

Semi - adjustable articulator with a calibration system.

Condylar Angle(0°~50°) / Bennett Angle (Fixed at 15°)



## CSA Series Articulators feature CALIBRATION System





Accurate & Precise results!



Lock





Magnetic Plate











Calibration System with adjustable condylar inclination and ergonomic design

CSA400 articulator is the perfect choice for andlysis, treatment planning, and completion of even the most complex occlusal treatment.

It combines the advantages of easy use, simplicity and ability to tilt the instrument back at a 30 degree angle.

## CSA 300 ARTICULATOR

Economical and Precision Easy use to tilt the instrument back at a 30 degree angle. Fully open access from the lingual position for easy access during denture setup so excellent to check lingual position when set array. Convenient attaching & removing cast with magnetic mounting plate.

### average movement articulator

Condylar Angle(fixed at 25°) / Benntt Angle(fixed at 15°)

## **CSA Series Articulators feature CALIBRATION System**





## CSA Accessory(Separately Sales)







CSA-36 Incisal Table (ADJUS TABLE)

CSA-2500 Curve of Spee

Average Mounting Jig (LOWER)

System

Precision

Operation

CSA300 articulator has been designed to be used accurately, safely and conveniently in all sectors of dental prosthetics such ad occlusion, treatment, analysis and so on.

back at a 30 degree.





Accurate & Precise results!







CSA-1600 Occlusal Plane Guide(UPPER)



CSA-1700 Cori Template

Average-valued type of Articulator Calibration system and Ergonomic desgin

It combines the advantages of easy use, simplicity and ability to tilt the instrument

## **FACE BOW**

### Product no. FB-1500

CORI FACE BOW is designed to figure out the occlusion relation with easy and accuracy.

One touch screw of Universal joint is perfect for transfer with simple & accuracy

## **ONE TOUCH SYSTEM**



Universal joint with Bite fork



Bow arm & Orbitale indicator



Transfer Jig



CSA600 With Face Bow



CSA400 With Face Bow



It is designed ergonomically so that the relationship of occlusion can be figured out easily.

The joints of Bite fork can be fixed by One touch screw at once and it can move left, right and centre freely.

Convenient operation by the simple designed structure.



### Please try CORIDENT SIMPLEX ARTICULATOR to all precision prosthodontics.

## Product no. E-200





## • Simple to use.



• Accessible to occlusion adjustment.

• Allowing vertical, hinge, lateral and slightly retrogressive movements.


Remounting metal magnetic caps.

• Economical For convenience, even if incisal pin is removed, the upper frame can be maintained in a level position by using the rear adjusting screw.

## SINGLE PIN

Designed dlaborately and has a unique shape in order to prevent distortion during pin(modeling) process.

SINGEL PIN DIMENSION(mm)				
L	L1	L2	D1	D2
21.00	16.65	5.00	1.80	2.00

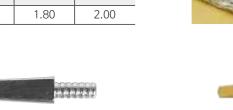



F-202 Material : Zinc nickel plated Packing: 1,000/box





F-203 Material : Zinc brass plated Packing: 1,000/box

## **TWIN PIN**


## Product no. F-210

Male and Female parts are perfectly fitte without any torsion. Male pin is combined very softly with Female(Sleeves). The main materials is Zinc and Nickel plated. Each rubber cap is packed individually for user's convenience.



Composition : 1,000 Sets/Box Male Pins 1.000 Pcs Female (Sleeves) 1.000 Pcs Rubbers 1,000 Pcs



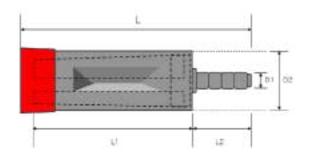

► F-210(TWIN PIN 16mm)			DIME	NSION(mm)
L	L1	L2	D1	D2
21.00	16.65	5.00	1.80	2.00

► F-210-12SHORT(TWIN PIN 12.5mm)	DIMENSION(mm)
----------------------------------	---------------

L	L1	L2	D1	D2
21.00	16.65	5.00	1.80	2.00

## **EURO PIN**

## Product no. F-300 Materials : Zinc, Nickel plated. Pin: 1,000pcs Plastic sleeve : 1,000pcs Rubber cab: 1,000pcs / Box 01

Introducing EURO PIN

- Well received in Europe.
- The female part has effective structure divided into two sleeves.
- EURO PIN has the function of compensating the hardening and expansion of anhydrite.

## **PLUS PIN**

### Product no. F-200

Materials: Zinc, Brass plated Pin : 1.000 EA Plastic sleeve : 1,000 EA Rubber cap: 800 EA / Box

## **EXCELLENT REMOVABLE PIN**

1. Plus pin wall perpendicularly (0 degree) parallel to female sleeve allows precision.

2. Miniature cross shape pin will prevent distortion.

3. Black colored rubber cap help easily find the position of pin on the bottom of poured stone cast.

4. Plus pin require only minimum space and it'll be easy to trim die and cast.




The pin-head part is more extended and has a thin round plate; therefore it supports the pin stably, distinguishing it from other products.

Thin structure of this part is for slight movement of it's sleeve with reflecting and compensating any change, caused by the expansion of the anhydrite.

► F-300 EURO PIN DIMENSION(m				NSION(mm)
L	L1	L2	D1	D2
20	13.7	5.3	2	5.9



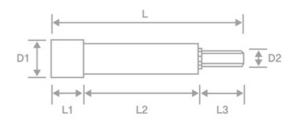

▶ PLUS PIN				DIMEN	SION(mm)
L	L1	L2	L3	D1	D2
21.00	3.20	13.00	5.00	5.00	2.0

## TELESCOP SOFT PIN


### Product No. F-110 & F-110-13

Materials: Zinc. Nickel plated Pin: 1,000 EA Plastic Sleeve: 1,000 EA/ Box

- Female sleeve of Telescop Soft pin is made of soft plastics.
- This soft female part provides users feeling tight in working with pin.
- The united shape of sleeve-body and bottom-cap reduces users' bother to put the cap on the pin.
- The specially thin parts in the united bottom cap helps to be separated easily from sleeve body.
- Female's keeping hole can catch male securely.

### ► F110 TELESCOP SOFT PIN

L	L1	L2	D1	D2
21.5	15.4	5	2	5.34

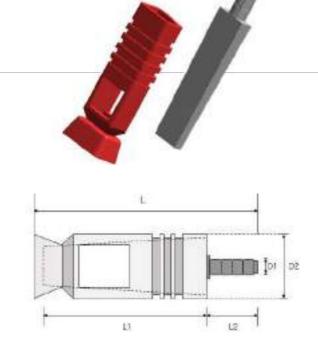

### ► F110-13 SHORT TELESCOP SOFT PIN

			-	
L	L1	L2	D1	D2
18.9	12.95	5	2	5.15

## **IMPLANT TRAY**

## Product No. B200 & B200-01

Metal Tray: 4 Pcs. / Set LARGE TRAY-UPPER & LOWER Parts MEDIUM TRAY-UPPER & LOWER Parts Plastic Cover : 210 EA/Set 105 ea. For L1~L7 parts-15 ea. of each parts 105 ea. For M1~M7 parts-15 ea. of each parts



- CORIIMPLANT TRAY is harmless to human with firm aluminum body which is strong, secure and can be disinfected in the high temperature.
- CORIMPLANT TRAY's transparent plastic cover are composed of 7 pieces which allow to see the screw zone and while get impression and to make the hole in the exact position of the cover. Also the each cover parts can be separated from each other and tray for user's convenience.
- It is convenient to get the impression without making individual try since the plastic part can be detached from the implant screw.
- It is more economic than individual tray and disposable tray.







## **PIN SYSTEM**

## Product No. F-700



## | PIN SYSTEM + Mounting ) Only Takes "30 minutes"

- Pin work and Mounting can be proceeded at the same time Existing Pin System takes 8 hours to finish, but CORI Pin system only takes 30minutes.
- It is compatible with other Articulators.
- Time saving, higher accuracy, and more convenient.
- Pin mold and articulator can be separated, so it is easy to install on surveyors.









